



Breffni Mobility
Ballydurrow, Virginia, Co Cavan

Telephone: 0871 22 00 00
 0873 29 08 40
 0498 54 84 54
fax: 049 854 8454

Dublin: 087 121 50 00
Galway: 087 121 60 00

www.breffnimobility.com
E mail: info@breffnimobility.com



Otolift Ruby Stairlift

The compact and
user-friendly stairlift




The Ruby Stairlift

The Ruby stairlift from Otolift is the most compact stairlift on the market. The Ruby is a typical example of Dutch design: a perfect combination of astute ergonomics, user-friendliness and neat design. The slim shape makes the stairs optimally accessible to other members of the household.

The Ruby stairlift is distinguished by its functional design. It guarantees comfort and safety for you and the other members of your household. The Ruby stairlift is mounted along the wall, without narrowing the treads. This is partly thanks to the uniquely compact track system, which consists of two slim tubes and is comparable to a standard handrail as regards shape and size. In most situations, the upper tube can also be used as a handrail. A great deal of attention has been paid to the styling of this stairlift. The result is an especially subtle, timeless and neutral design that fits in with any interior.

The Ruby stairlift
gives you renewed
freedom



Discovering that you can no longer do certain things that you used to do without thinking is concerning. You have to learn to cope with your limitations and accept that you can no longer do everything under your own steam. You have to think of possibilities and solutions. With the Ruby stairlift from Otolift, you will be able to continue living in your familiar surroundings and regain part of your lost freedom and independence.

STABLE, COMFORTABLE RIDE QUALITY


The unique track system ensures a stable, quiet ride quality. The double track guarantees 100% stability – the stairlift simply cannot topple. The Ruby provides the user with pleasant and stable transport between floors, with no rocking or vibration. A 'soft' start and stop contributes to the comfort of the ride – the lift starts and stops shock-free. Thanks to the continuous transmission between the various speeds, the ride quality of the ascent and descent is smooth.

SAFE AND RELIABLE

Safety, comfort and reliability are the foremost aspects of this design. The comfortable seat and back support ensure that you feel secure and at ease. The Ruby is provided with many safety features, such as a safety belt. If the lift comes up against the slightest resistance, it stops immediately. And on arrival at your destination, you simply give the seat a quarter turn, thus forming a barrier across the stairwell and allowing you to step off safely and easily.



↑ The least resistance will stop the stairlift immediately.

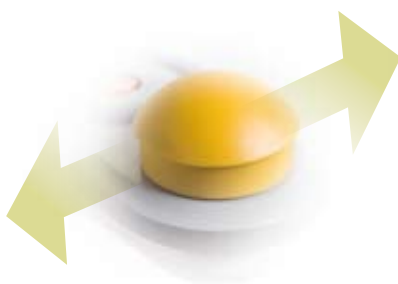

↑ The comfortable chair ensures that you move safely and relaxed between floors. You set the stairlift in motion using the yellow ball on top of the armrest.



↑ Getting on and off at the top of the stairs is done safely by simply giving the seat a quarter turn.







↑ Thanks to its subtle design, the Ruby stairlift fits perfectly in with the interior. And the slim track means that the stairs remain freely accessible to visitors and other members of the household. The double-track system makes the lift 100% stable, so that the lift chair cannot topple.



↑ Folded away, the stairlift is very compact. The mechanism is cleverly concealed.



↑ The running track is similar to a handrail in size.

USER-FRIENDLY

Operation of the Ruby stairlift excels in its simplicity. The compact chair provides optimum convenience: all the elements, such as the footrest, armrest and seat, fold easily away. When not in use, the stairlift takes up a minimum of space. You set the lift in motion using the yellow operating ball on top of the armrest. The unique ball shape enables anyone to operate the lift easily and without effort. The Ruby is supplied with remote control as standard for summoning and sending the lift, with no cables necessary.

RELIABLE IN A POWER CUT

The Ruby chairlift is powered by two batteries that are automatically charged when the lift is parked on its usual spot. The charging is done automatically via the two charging stations positioned on the track. You are not dependent on the electricity supply. In the event of a power cut, the stairlift will function undisturbed. So you need never fear that you will get stuck halfway – you will still be able to go up or down several times if the electricity supply is disrupted for whatever reason.

SIMPLE TO INSTALL

Installation takes place quickly and simply. The lift is securely fastened on to your stairway, without the need to drill into the wall. The Ruby can be installed on practically any type of stairway, and is made to measure. Otto Ooms offers a range



of optional extras for adapting the Ruby stairlift to your personal needs or your specific living circumstances. For the lift chair upholstery, there is a choice of twelve colours and two types of covering; there is a choice of two colours of paintwork for the lift. Custom solutions, such as telescope or foldaway extension and a parking function, make the Ruby stairlift ideally suitable in restricted spaces.



← There are various means available for keeping the passageway free at the foot of the stairs



The Ruby Stairlift: refined in the details.

TEN BENEFITS OF THE RUBY STAIRLIFT:


- The stairway remains freely accessible
- 100% stable
- Subtle style
- Safe and comfortable
- Reliable
- Wireless remote control
- Soft start and stop
- Will operate during a power cut
- Simple to install on any type of stairway
- Very quiet

TWELVE CHAIR COVERINGS

The coverings of the seat and backrest are available in eight materials and four imitation leather choices.



Violet

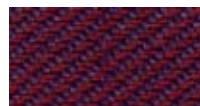

Burgundy



Navy blue



Grey



Brown



Green



Red



Light brown



Hazelnut

Imitation leather



Burgundy

Imitation leather



Black

Imitation leather




Midnight blue

Imitation leather

TWO COLOURS OF PAINTWORK



RAL 7035



RAL 9001



TECHNICAL SPECIFICATIONS

Lifting capacity 130 kg

Speed Variable to max 0.15 m/sec with soft start and stop

Motor 300 Watt

Electricity supply 230 Volt 1 phase socket in the vicinity of the stair for the battery charger

Controls On the chair, in the left or right armrest. Call and send switch at the top and bottom of the stairs.

Rail Running track supplied with powder coating, mounted on the stair treads

Special Design

- Model suitable for children
- Model suitable for stiff limbs / joints
- Partly weather-resistant model

Many other applications possible upon request

Safety Pressure-sensitive features under the footrest and chair, the side panels, and driving mechanism. Rollup safety belt

Extras Key switches and various personalised adjustments

GUARANTEE AND SERVICE

Otolift stairlifts are innovative quality products with a long lifecycle. Servicing and any repairs are carried out quickly and expertly by our own maintenance team. Our service department is open all year round, 7 days a week.

